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
AP721792-10_Thermometer Bracelet Leroux User manual (Eng)

It is recommended to read this user manual carefully and keep it for future reference.

Product overview:

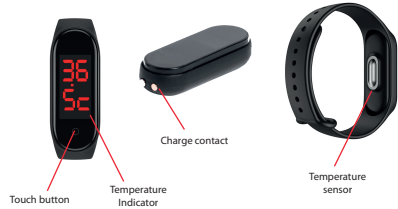Specifications:

- Charging Voltage: 5V
- LED Screen: 0.96 inches
- Li-polymer battery 55 mAh, 3.7V
- Watch working time: 30 days.
- Degree of Protection: IP54
- Temperature range: 34-42.5°C

Charging:

Use the supplied cable for charging. Connect the USB cable to your computer / AC adapter and connect the other end of the cable to the charge contact of the device.

Watch operation:

Long press the touch button for 5 seconds to enter setting mode. The screen will show hour digits. Press the touch button to adjust.

Long press the button for 5 seconds to set the minutes and press to adjust.

When the correct hour and minutes appear on the screen, hold down the button to save the time.

Thermometer operation:

Press the touch button until the screen shows the temperature.

The letter "C" situated on the bottom right corner will start to flash. This indicates that the body temperature reading is in process.

When the temperature is reached, the letter "C" will stop flashing and remain on the screen.

Long press the touch button to show the temperature in °F.

Note: When you put the bracelet on, the sensor will take a few minutes to read your body temperature. When the body temperature exceeds 37.5°C a flashing red light will appear in the top left corner. This warning light will flash for 10 seconds.

Note: If the temperature is lower or higher than the watch temperature range (34-42°C) LO or HI respectively will appear on the screen.

Warnings:

1. Do not expose this product to rain or damp conditions.
2. Prevent this item from falling as it could cause serious damage.
3. Do not dismantle, repair or modify this product by yourself, contact an authorized specialist.
4. In order to avoid any risk, do not dispose of this product inappropriately or throw it into the fire as there is a built-in Lithium battery.
5. Clean it with a dry cloth.
6. Do not place this item near sources of heat or expose to direct sunlight.
7. Keep out of reach of children.
8. Do not insert any metal objects into this device as there is a risk of causing a short circuit in the unit.

This symbol both in the product and its user manual indicates that at the end of the device's electrical life, it should be recycled separately from your household waste. There are adequate means to dispose of this material for proper recycling. For more information, contact your local authorities.

Hereby, Miko Catal Importaciones S.L. declares that this radio equipment is in compliance with Directive 2014/53/EU. The full text of the UE declaration of conformity is available at: www.mikoto.eu (search on reference).

Distributor: Cool Collection, 1476, Budapest, POB:176, Hungary
Registration nr.: 01-09-276783, Made in China

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